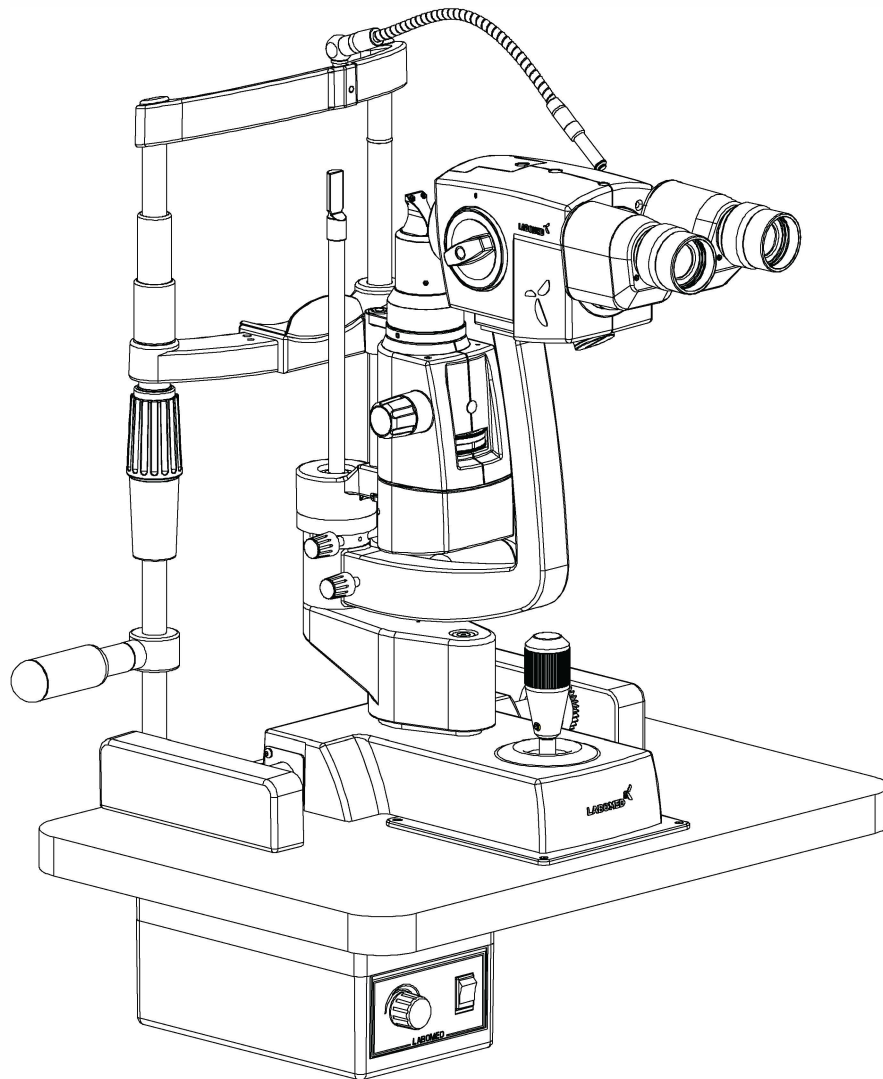


eVO 300

Slit Lamp

User Manual



eVO 300 is a registered Trade Mark of Labotech Microscopes India Pvt. Ltd.

Labomed is the registered Trade Mark of Labo America Inc.

All other trade marks are property of their respective owners.

The information contained in this document was accurate at the time of publication. Specification subject to change without notice.

Labotech and Labo America Inc. reserves the right to make change in the product described in this manual without notice and without incorporating those changes in any products already sold.

ISO 13485 Certified - LABOMED products are designed and manufactured under quality management system conforming to ISO 13485 requirements.

No parts of this publication may be reproduced, stored in the retrieval system, or transmitted in any form or by any means, electronic, mechanical, recording, or otherwise without the prior written permission of LABOTECH.

Caution: Federal law restrict this devices to sale by or on the order of a licensed physician. Rx only.

TABLE OF CONTENTS

Warnings and Cautions	4
Symbol Information	6
Introduction	7
Indications for Use	7
Contraindications	7
Setup	8
Parts Identification	8
eVO 300 Package Contents	8
Accessories	8
Unpacking & Installation	9
Application of Input Power	11
Disconnection of Input Power	11
Setup Of Camera	12
Instructions for Use	13
Operation	13
Slit Length	14
Filter Dial	14
Slit Rotation	14
Cleaning & Maintenance	15
External Cleaning	15
Forehead / Chin rest Preparation	15
Cleaning the Glide Plate	15
Changing the Halogen Bulb	15
Changing the LED	16
Fuse Replacement	17
Replacement of Chin Rest Light Indicator	17
Troubleshooting	18
System Configuration	19
Chart of Common Errors	20
Specifications	20
Disposal	21
Software Revision	21
Guidance Tables	22
Warranty	26

WARNINGS AND CAUTIONS

LABOMED is not responsible for the safety and reliability of this instrument when:

- Assembly, disassembly, repair, or modification is made by unauthorized dealers or persons.
- The instrument is not used in accordance with this user manual.

A **WARNING** is an instruction that draws attention to the risk of injury or death.



WARNING: UNITED STATES FEDERAL LAW AND EUROPEAN REGULATIONS REQUIRE THAT THIS DEVICE BE PURCHASED ONLY BY A PHYSICIAN OR A PERSON ACTING ON BEHALF OF A PHYSICIAN .

WARNING: THIS INSTRUMENT SHOULD BE USED IN STRICT ACCORDANCE WITH THE INSTRUCTIONS OUTLINES IN THIS USER'S GUIDE. THE SAFETY OF THE OPERATOR AND THE PERFORMANCE OF THE INSTRUMENT CANNOT BE GUARANTEED IF USED IN A MANNER NOT SPECIFIED BY LABOMED TECHNOLOGIES.

WARNING: DO NOT REPAIR OR SERVICE THIS INSTRUMENT WITHOUT AUTHORIZATION FROM THE MANUFACTURER. ANY REPAIR OR SERVICE TO THIS INSTRUMENT MUST BE PERFORMED BY EXPERIENCED PERSONAL OR DEALERS WHO ARE TRAINED BY LABOMED OR SERIOUS INJURY TO THE OPERATOR OR PATIENT MAY OCCUR.

WARNING: MODIFICATIONS TO THIS INSTRUMENT ARE NOT ALLOWED. ANY MODIFICATION TO THIS UNIT BE AUTHORIZED BY LABOMED OR SERIOUS INJURY TO THE OPERATOR OR PATIENT MAY OCCUR.

WARNING: IF THIS INSTRUMENT IS MODIFIED, APPROPRIATE INSPECTION AND TESTING MUST BE CONDUCTED TO ENSURE CONTINUED SAFE USE OF THIS INSTRUMENT.

WARNING: TO AVOID RISK OF ELECTRIC SHOCK, THIS EQUIPMENT MUST ONLY BE CONNECTED TO A SUPPLY MAIN WITH PROTECTIVE EARTH OR DAMAGE TO THIS INSTRUMENT AND/OR INJURY TO THE OPERATOR OR PATIENT MAY OCCUR.

WARNING: ENSURE THAT THE VOLTAGE APPLIED TO THE UNIT IS THE SAME AS THE VOLTAGE THAT IS INDICATED ON THE DATA PLATE OR DAMAGE TO THE UNIT MAY OCCUR.

WARNING: THIS INSTRUMENT MUST BE PLUGGED INTO AN OUTLET WITH AN EARTH GROUND. DO NOT REMOVE OR DEFEAT THE EARTH GROUND CONNECTION ON POWER INPUT CONNECTOR OR THE UNIT'S POWER CORD OF THIS INSTRUMENT OR DAMAGE TO IT AND/OR INJURY TO THE OPERATOR OR PATIENT MAY OCCUR.

WARNING: THE EQUIPMENT OR SYSTEM SHOULD NOT BE USED ADJACENT TO OR STACKED WITH OTHER EQUIPMENT AND THAT IF ADJACENT OR STACKED USE IS NECESSARY, THE EQUIPMENT OR SYSTEM SHOULD BE OBSERVED TO VERIFY NORMAL OPERATION IN THE CONFIGURATION IN WHICH IT WILL BE USED.


WARNING: THIS INSTRUMENT IS NOT SUITABLE FOR USE IN THE PRESENCE OF FLAMMABLE ANESTHETIC MIXTURES, SUCH AS OXYGEN OR NITROUS OXIDE.

WARNING: BECAUSE PROLONGED INTENSE LIGHT EXPOSURE CAN DAMAGE THE RETINA, THE USE OF THE DEVICE FOR OCULARE EXAMINATION SHOULD NOT BE UNNECESSARILY PROLONGED, AND THE BRIGHTNESS SETTING SHOULD NOT EXCEED WHAT IS NEEDED TO PROVIDE CLEAR VISUALIZATION OF THE TARGET STRUCTURES. THIS DEVICE PROVIDED WITH FILTERS THAT ELIMINATE UV RADIATION (<400NM) AND, WHENEVER POSSIBLE, FILTERS SHORT- WAVELENGTH BLUE LIGHT (<420NM).

WARNING: THE USE OF ACCESSORIES OR CABLES OTHER THAN THOSE SPECIFIED, WITH THE EXCEPTION OF THOSE SOLD BY THE MANUFACTURER AS REPLACEMENT PARTS FOR THE INTERNAL COMPONENTS, MAY RESULT IN INCREASED EMISSIONS OR DECREASED IMMUNITY OF THE EQUIPMENT OR SYSTEM.

WARNINGS AND CAUTIONS (continued)

A CAUTION is an instruction that draws attention to the risk of damage to the product.

 **CAUTION:** THE INTERNAL CIRCUITRY OF THE INSTRUMENT CONTAIN ELECTROSTATIC SENSITIVE DEVICES (ESD) THAT MAY BE SENSITIVE TO STATIC CHARGES PRODUCED BY THE HUMAN BODY. DO NOT REMOVE THE COVERS WITHOUT TAKING PROPER ESD PRECAUTIONS.

CAUTION: DO NOT USE SOLVENTS OR STRONG CLEANING SOLUTIONS ON ANY PART OF THIS INSTRUMENT, AS DAMAGE TO THE UNIT MAY OCCUR SEE THE CARE AND MAINTENANCE SECTION FOR DETAILED CLEANING INSTRUCTIONS.

CAUTION: MEDICAL ELECTRONIC EQUIPMENT NEEDS SPECIAL PRECAUTIONS WITH RESPECT TO ELECTROMAGNETIC CHARGE (EMC) AND NEEDS TO BE INSTALLED AND SERVICED ACCORDING TO THE EMC INFORMATION PROVIDED IN THE ACCOMPANYING DOCUMENTS.

CAUTION: PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT CAN AFFECT MEDICAL ELECTRICAL EQUIPMENT.

CAUTION: THIS INSTRUMENT IS NOT TO BE USED NEAR HIGH-FREQUENCY EMITTING SURGICAL EQUIPMENT.

CAUTION: THIS INSTRUMENT IS NOT INTENDED TO BE CONNECTED TO EQUIPMENT OUTSIDE THE CONTROL OF LABOMED OR MUST BE TESTED TO AN APPLICABLE IEC OR ISO STANDARDS.

CAUTION: CHIN REST PAPER-FOR SINGLE USE ONLY. REPLACE WITH EACH PATIENT.

EXPLANATION OF SYMBOLS

The following symbols appear on the instrument:



Caution symbol indicating important operating and maintenance Instructions that are included in this User's Guide.



Type B Applied Part



Alternating Current Power



Protective Earth



Connection ON/OFF



Manufacture

REF

Catalog Number

S/N

Serial Number



Waste of Electrical and Electronic Equipment



Compliance to medical Device Regulation EU MDR 2017/745



Accompanying Document must be consulted.

Authorized Representative in European Community.



Fragile Contents in Shipping Container- Handle with care



Keep Dry- Package shall be kept away from rain.



This way Up- Indicates correct upright position of package.



U.S. Federal law restricts this device to sale by or on the order of a physician.



Medical Device


INTRODUCTION

Congratulations on your purchase of the eVO 300 Slit Lamp.

This User Manual is designed as a training and reference manual for the operation and maintenance of the instrument. We recommend that you read it carefully prior to use and follow the instructions to ensure optimum performance of your new instrument. Properly trained eye care professionals such as ophthalmologists, optometrists, opticians and eye care technicians should operate this instrument.

Please retain this manual for future reference and to share with other users. Additional copies can be obtained from your authorized LABOMED dealer or from the LABOMED Customer Service Department at:

Tel: (510) 445-1257
Fax: (510) 991-9862
Email: sales@laboamerica.com

 Any serious incident associated with this device should be reported to Labo America, Inc. USA and to your national health agency.

INDICATIONS FOR USE

The eVO 300 Slit Lamp is an AC-powered slit lamp biomicroscope that is intended for use in examining the anterior segment, from the corneal epithelium to the posterior capsule. It is used to aid in the diagnosis of diseases or trauma, which affect the structural properties of the anterior segment of the eye.

CONTRAINDICATIONS

None.

FEATURES AND FUNCTIONS

PART IDENTIFICATION

1. On/Off Switch
2. Illumination Level Control
3. Glide Plate
4. Joystick for horizontal and vertical movement
5. Bulb Access Door
6. Filter Dial
7. Slit Length Dial
8. Breath Shield Mount
9. Eyepieces
10. Focusing Ring
11. Fixation Light
12. Magnification Dial
13. Microscope Lock Knob
14. Slit Rotation Scale
15. Slit Width/Rotation Knob
16. Illumination Arm Lock Knob
17. Microscope Arm Lock Knob
18. Instrument Base Lock Knob
19. Guide Rail Cover
20. Geared Rollers
21. Inbuilt tube for CCD Camera Attachment
22. Chin Rest Assembly
23. Power Supply Assembly
24. Table Top
25. Patient Handles
26. Yellow Filter Knob

eVO 300 PACKAGE CONTENTS

1. eVO 300 Slit Lamp HL (8126200)
2. eVO 300 Slit Lamp LED (8126201)
3. User's Guide (8126200-795)

ACCESSORIES

1. Focusing Rod (8124300-801)
2. Hex Wrench - 3mm (LK-003)
3. Hex Wrench - 4mm (LK-004)
4. Hex Wrench - 5mm (LK-005)
5. Dust Cover (PKG-107)
6. Halogen Lamp, Main (8124300-401)
7. Rack Covers (8124100-206)
8. Replacement Fuses (8124900-900)
9. Chin Rest Paper - 1 pack (8126000-786)



SETUP

UNPACKING AND INSTALLATION

1. Open the outside shipping box and remove the four (3) inner boxes.
2. Remove the User's Guide and follow it for Slit lamp installation.
3. Open the box with the Table Top and Electronics in it. Refer to Figure-1
4. Remove the Table Top from the box and install the Table Top onto the instrument stand or table, and secure it in. Place as indicated in the user guide for the stand.
5. Open the box with the chin rest assembly and remove it. Refer to Figure-2.
6. Using the 3mm Hex Wrench, connect the Ground Wire from the Power Supply to the Chin Rest Assembly using the Ground Screw provided (torque to 7.8 N.m). Refer to Figure-3.
7. Using the 5mm Hex Wrench, remove the two Screw from the button of the Table Top and attach the Chin Rest Assembly to the Table Top using these screws. Refer to Figure-4.



Figure-1 Table Top and Electronics



Figure-2 Chin Rest Assembly



Figure-3 Chin Rest Ground

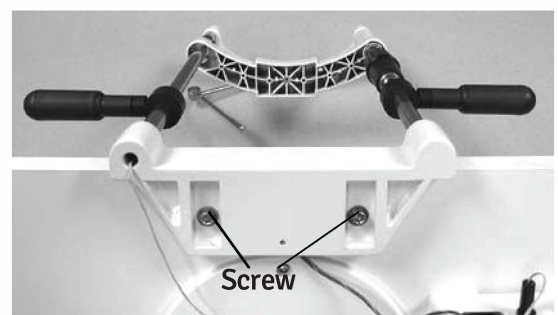


Figure-4 Chin Rest Attaching Screw

SETUP (continued)

UNPACKING AND INSTALLATION (continued)

8. Attach the Fixation Light Wire from the Chin Rest Assembly into the connector on the back of the Power Supply Assembly. Refer to Figure-5.

9. Open the box with the Microscope Assembly. Remove the Microscope Assembly, Slit lamp Assembly, and Accessories. Refer to Figure-6.

10. Install the Slit Lamp Assembly onto the tracks of the Table Top and slide the Guide Rail Covers around the tracks. Refer to Figure-7.

11. Attach the Base Lamp Wire to the back of the Power Supply Assembly. Refer to Figure-5.

12. Using the 4mm Hex Wrench, adjust the Patient Handles by loosening the Allen Cap screws that are securing them to the Chinrest Posts. Slide the Patient Handles up or down to the desired height, and secure them in place by tightening the Allen Cap Screws. Refer to Figure-8.

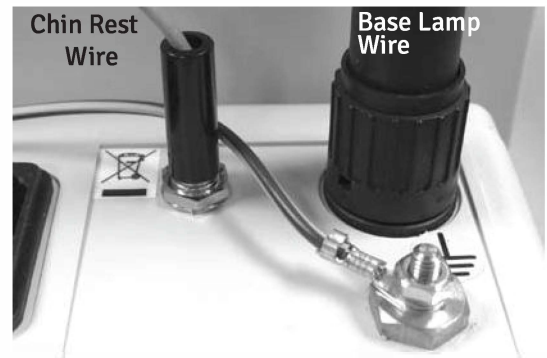


Figure-5 Connections

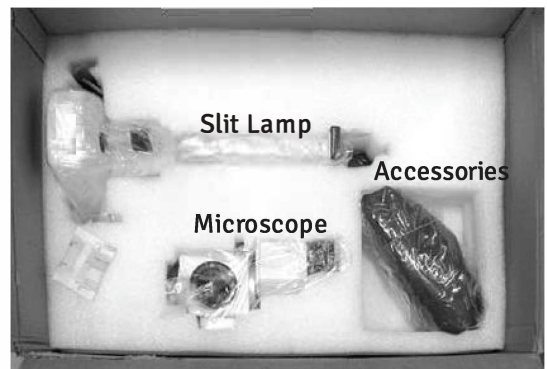


Figure-6 Microscope, Base, Accessories



Figure-7 Install Base

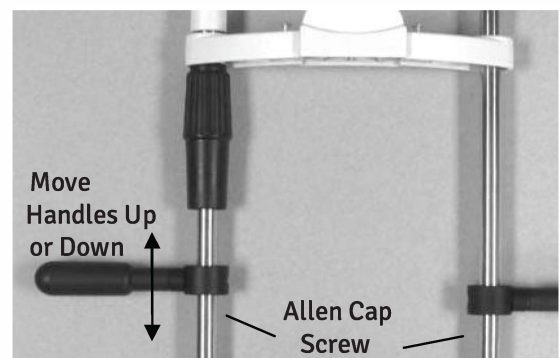


Figure-8 Adjusting Patient Handle Height

SETUP (continued)

UNPACKING AND INSTALLATION (continued)

13. Install the Microscope Assembly onto the top of the Arm by sliding it into position, making sure it is up against the stop. Then, tighten the Lock Knob located on the right side of the Microscope Assembly. Refer to Figure-9.

Note: Do not adjust the microscope stop knob behind the base of the microscope, or the vertex distance will cause misalignment of focus and require re-calibration of the slit lamp assembly.

14. Remove the accessories and store them in an appropriate place so that when they are needed they will be available. Refer to Figure-10.

APPLICATION OF INPUT POWER

WARNING: CARE MUST BE TAKEN TO ARRANGE THE CABLES FOR THE ACCESSORIES SUCH THAT THEY DO NOT PRESENT A TRIPPING HAZARD TO THE EXAMINER OR A DANGER TO THE PATIENT.

WARNING: POSITION THIS INSTRUMENT SO THAT IT IS NOT DIFFICULT TO OPERATE THE DISCONNECTION DEVICE(PLUG).

1. After the unit is in its secure location, apply the correct input voltage to the instrument using the Power Cord from the Accessory Tray.

NOTE: The power inlet is located on the backside of the Power Supply Assembly.

2. Press down on the "I" located on the ON/OFF Switch. Refer to figure-11.

NOTE: The ON/OFF Switch will illuminate green when there is power to the unit. When the ON/OFF Switch is set to off, the green light will turn off.

DISCONNECTION OF INPUT POWER

1. At any time, the power switch can be set to OFF. The unit does not have a power down sequence. To terminate operation of this instrument, press the ON/OFF switch to the OFF position (O).

2. If this instrument is intended to be OFF for an extended period of time, it can be disconnected from power by detaching the power cord from the receptacle.



Figure -9 Install Microscope



Figure -10 Accessories

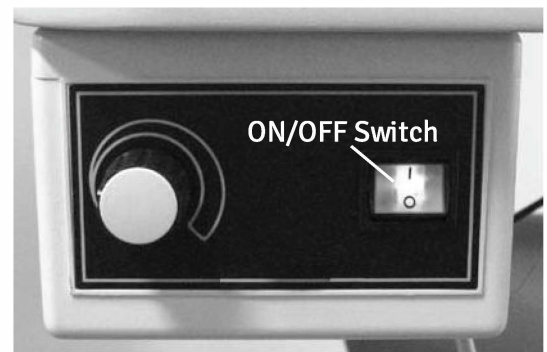


Figure-11 Power Supply Assembly

SETUP (continued)

UNPACKING AND INSTALLATION (continued)

Set Up of Camera (CCD / HD CCD Camera)

1. Remove the Check Nut from the Inlet Tube for CCD Camera shown as (A)., Refer Figure. 12.
2. Fix the Lens Adapter by rotating it Clock wise to the Inlet Tube for Camera shown as (B)., Refer Figure. 13.
3. Remove the Dust Cover from the Camera and fix the Camera by rotating it clockwise into the Lens Adapter. Shown as (C)., Refer Figure. 14.
4. Plug in the video cable into the Analog video Connector. Shown as (E)., Refer Figure. 15.
5. Plug in the Power Connector into the Camera Power inlet. Shown as (F)., Refer Figure. 15.
6. Connect the Video Cable into The TV/LED Video in Port to Live streaming.

Note: To Align the Image on Camera with Image in Eye Piece.

1. Loosened the Groove Allen Screw by using 1.6mm Allen Key. Shown as (D). Refer image 14.
2. Rotate the Camera to align the image on camera, once set then tighten the Groove Allen Screw.

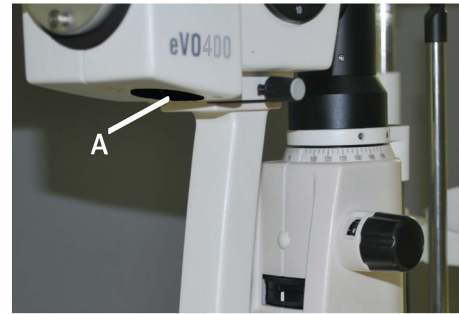


Figure 12

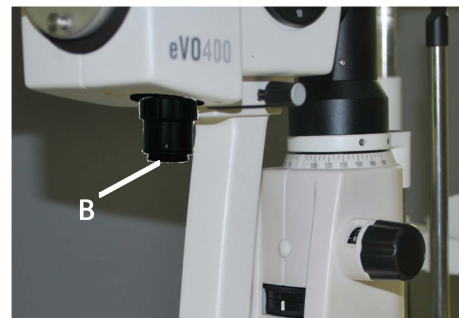


Figure 13

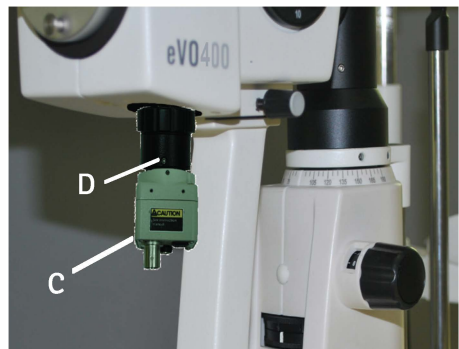


Figure 14



Figure 15

INSTRUCTIONS FOR USE

OPERATION

1. Turn on the power using the On/Off switch located on the front of the Power Supply Assembly. Brightness can be adjusted by rotating the illumination level knob.

NOTE: The maximum position is for intermittent use only. Continuous use will shorten the lamp/ LED life.

2. Insert the Focusing Rod in the pivot post of the instrument body to make rough IPD and focus adjustments.

3. Position the light onto the flat surface of the Focusing Rod and adjust the pupillary distance and focus of the eyepieces to suit the needs of the operator. Refer to Figure-12.

4. Using the Slit Width Knobs, adjust the projected slit so that the thinnest slit is shown on the Focusing Rod. Refer to Figures 12 and 15.

NOTE: The thinnest line will allow for greater accuracy.

5. Remove the Focusing Rod.

6. To position a patient, adjust the chinrest height by turning the Chinrest Elevation Handle on the post of the Chin Rest Assembly until the patient's canthus is in line with the canthus mark on the chin rest post. Refer to Figure-13.

7. Microscope elevation is adjusted by rotating the Joystick and observing the slit image through the Microscope Assembly until the slit is centered on the patient's cornea. Refer to Figure-14.

8. Move the slit lamp with the Joystick held firmly and slightly angled toward the patient, until the slit appears sharply on the cornea.

NOTE: The accuracy of this rough adjustment should be checked by the naked eye. The fine adjustment is performed while observing the slit through the microscope.

9. Tilt the Joystick, which is now held lightly at its upper end, until the slit appears sharply at the depth of the eye which is to be observed.

10. The horizontal motion of the base can be locked by tightening the Base Locking Screw. Refer to Figure-14.

NOTE: Lock the base whenever the lamp is not in use.

11. The slit width can be adjusted by rotating the Slit Width/Rotation Knob on either side of the instrument. Refer to Figure-15.

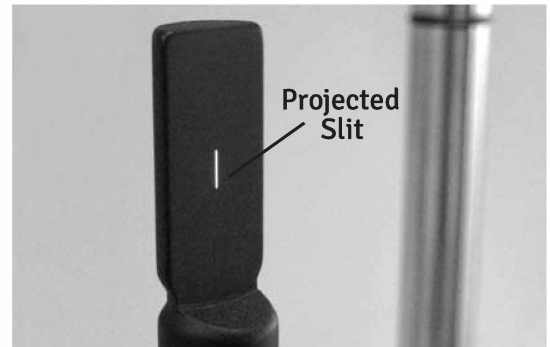


Figure-12 Focus on Slit



Figure-13 Adjust Patient Height

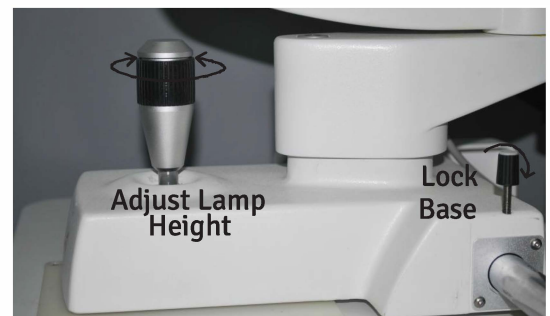


Figure-14 Adjust Height

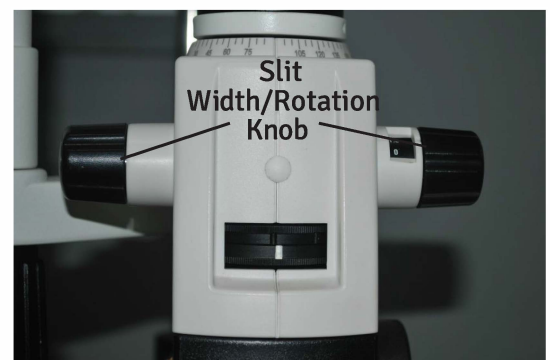


Figure-15 Adjust Slit Width

INSTRUCTIONS FOR USE (continued)

OPERATION (continued)

- 12. The angle between the illumination system and the microscope can be varied between 0° and 90° to either the left or to the right. Refer to Figure-16.
- 13. The illumination angle is indicated on the Scale of the slit lamp arm. Refer to Figure-17.
- 14. Magnification is altered by rotating the Magnification Dial on the Microscope Assembly. Refer to Figure-16.



Figure-16 Illumination Angle

SLIT LENGTH

The slit length is adjusted by rotating the Slit Length Dial. The dial has five stops for adjustments. They are 0.6, 5.8, 9, 13.5 mm diameter and continuous 1.5 to 12.0mm. They index into place. Refer to Figure-18.



Figure-17 Illumination Angle Scale

FILTER DIAL

The Filter Dial has four positions that index into place, and are color coded to indicate the active filter. Refer to Figure-18. The color coded index stops are as follows:

- Blue dot = Cobalt Blue
- Red dot = Heat Absorbing
- White dot = Open
- Green dot = Red-free



Figure-18 Filters and Slit Length

SLIT ROTATION

Slit rotation is achieved by grasping the Slit Width/Rotation Knob and twisting the slit body to the left or right. The degree of rotation is indicated by the Slit Rotation Scale above the slit body. Refer to Figure-19.

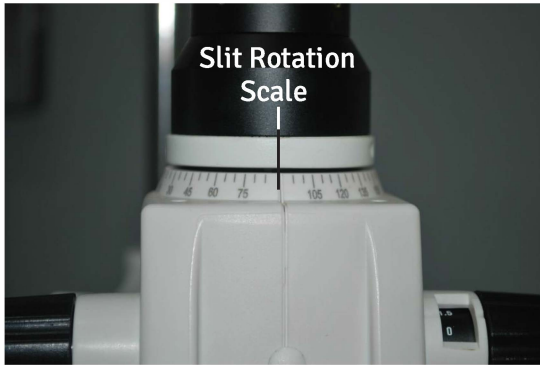


Figure-19 Slit Rotation Scale

CLEANING & MAINTENANCE

WARNING: RISK OF ELECTRIC SHOCK. ALWAYS DISCONNECT THE POWER CORD FROM THE WALL AND THE INSTRUMENT BEFORE PERFORMING ANY OF THE FOLLOWING CARE AND MAINTENANCE PROCEDURES.

EXTERNAL CLEANING

Clean the external surfaces of this instrument using a Clean, soft cloth moistened with a mild detergent solution (1 cc of liquid dish soap to one liter of clean, filtered water (filtered below 5 microns). Refer to Figure-20.



Figure-20 Cleaning Main Unit

FOREHEAD/CHIN REST PREPARATION

For hygienic reasons, wipe the forehead rest with an alcohol wipe and change the chin rest papers after each patient.

CLEANING THE GUIDE PLATE

If the Guide Plate is dirty it may cause a rough feeling when maneuvering the base of the slit lamp. Clean the Guide Plate with a soft cloth lightly dampened with a mild soap and water solution.



Figure-21 Bulb Holder

CHANGING THE HALOGEN BULB

WARNING: NEVER REMOVE A BULB THAT HAS RECENTLY BEEN IN USE AS IT WILL BE VERY HOT. WAIT UNTIL IT HAS COOLED AND USE GLOVES OR A THICK CLOTH WHEN HANDLING ANY HALOGEN BULB.

WARNING: NEVER TOUCH A HALOGEN BULB WITH BARE HANDS AS FINGERPRINTS WILL SHORTEN THE BULB LIFE.

1. Remove input power to the instrument.
2. Open the bulb door and pull out the bulb holder by pulling the metal pins on the holder. Refer Figure 21.
3. Swing the retaining clamp from the bulb by loosen the screw. Refer Figure 22.
4. Pull the bulb from the unit and replace the bulb with the correct bulb as indicated in the specifications section of the manual.

Note: Position the Bulb so the Cut Out in the metal collar of the Bulb lines up with the Notch in the lamp housing. Refer to Figure-24.

7. Move the Retaining Clamp back into its original position. Refer to Figure-22.
8. Place the bulb holder back on the position
9. Close the bulb door.

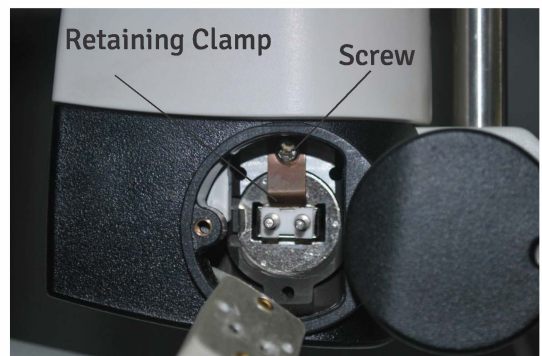


Figure-22 Retaining Clamp



Figure-23 Bulb

CHANGING THE LED

(For LED illumination only)

WARNING: NEVER REMOVE A LED THAT HAS RECENTLY BEEN IN USE AS IT MAY BE VERY HOT WAIT UNTIL IT COOLED.

WARNING: NEVER TOUCH A LED WITH BARE HANDS AS FINGERPRINT WILL SHORTEN THE LED LIFE.

1. Remove input power to the instrument.
2. Remove the Circlip by using circlip plier. shown as "A" Refer Figure 24.
3. Lift up the illumination housing carefully.
4. Remove the Lamp Cover by unscrewing the Four Screws by using screwdriver shown as "B". Refer Figure 25.
5. Unscrew the two M3 allen screws to disassemble the Heat Sink by using 2.5 Allen wrench shown as "C" Refer Image 26.
6. Unscrew the LED from Heat Sink by unscrewing the Two Screws shown as "D". Refer to Figure 27.
7. Replace New LED and follow Reverse to complete assembly.



Figure 24

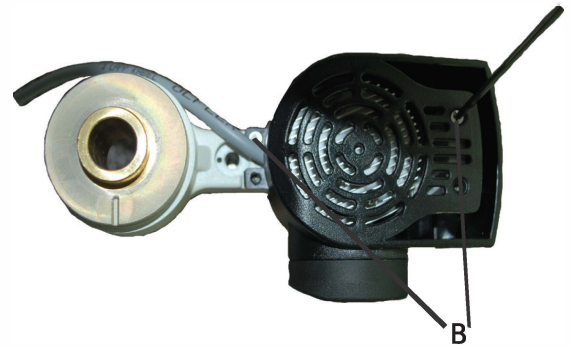


Figure 25



Figure 26

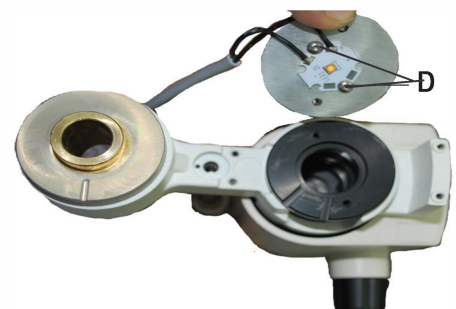


Figure 27

FUSE REPLACEMENT

Replace the fuses in the Power Input Module with the fuses indicated in the **Specifications** section of this manual.

1. Remove input power to the instrument.
2. Press down on the top tab in the middle of the Power Input Module to release the Fuse Holder, and gently pull out the Fuse Holder by gripping the two small tabs. Refer to Figures 25 and 26.
3. Open the Door to the Fuse Holder by pulling it down. Refer to Figure-25.

NOTE: The Fuses will pop up when the door is open, making removal easier.

4. Install new fuses into the Fuse Holder that is indicated in the Specification section of this manual.
5. Install the Fuse Holder by closing the door, and pushing the Fuse Holder back until it snaps into place.

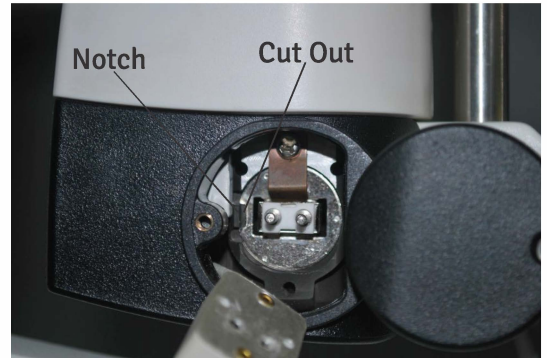


Figure-24 Notches

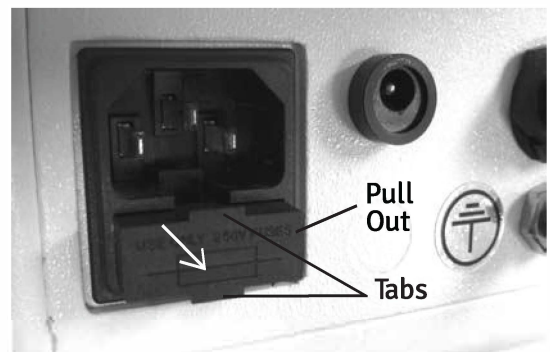


Figure-25 Pull Out

Replacement of Chin Rest Light Indicator

1. Unscrew the LED cap refer fig. 27.
2. Pull out the defective LED and replace with new LED. refer fig.27a. Thread in the LED cap back.

NOTE: If New LED does not light up after switching on. remove it and re-fix after changing its polarity.

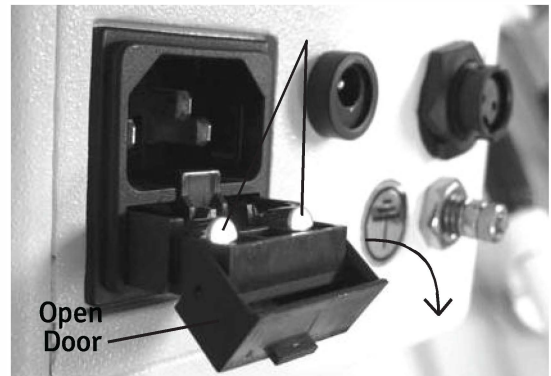


Figure-26 Open Fuse Door



Figure 27 Changing Target Light LED



Figure 27a Changing Target Light LED

TROUBLESHOOTING

The following chart outlines some common issues with the eVO 300 Slit Lamp and some steps you can take to correct the issue. If problems persist, please contact the LABOMED as listed in the Introduction section of this manual.

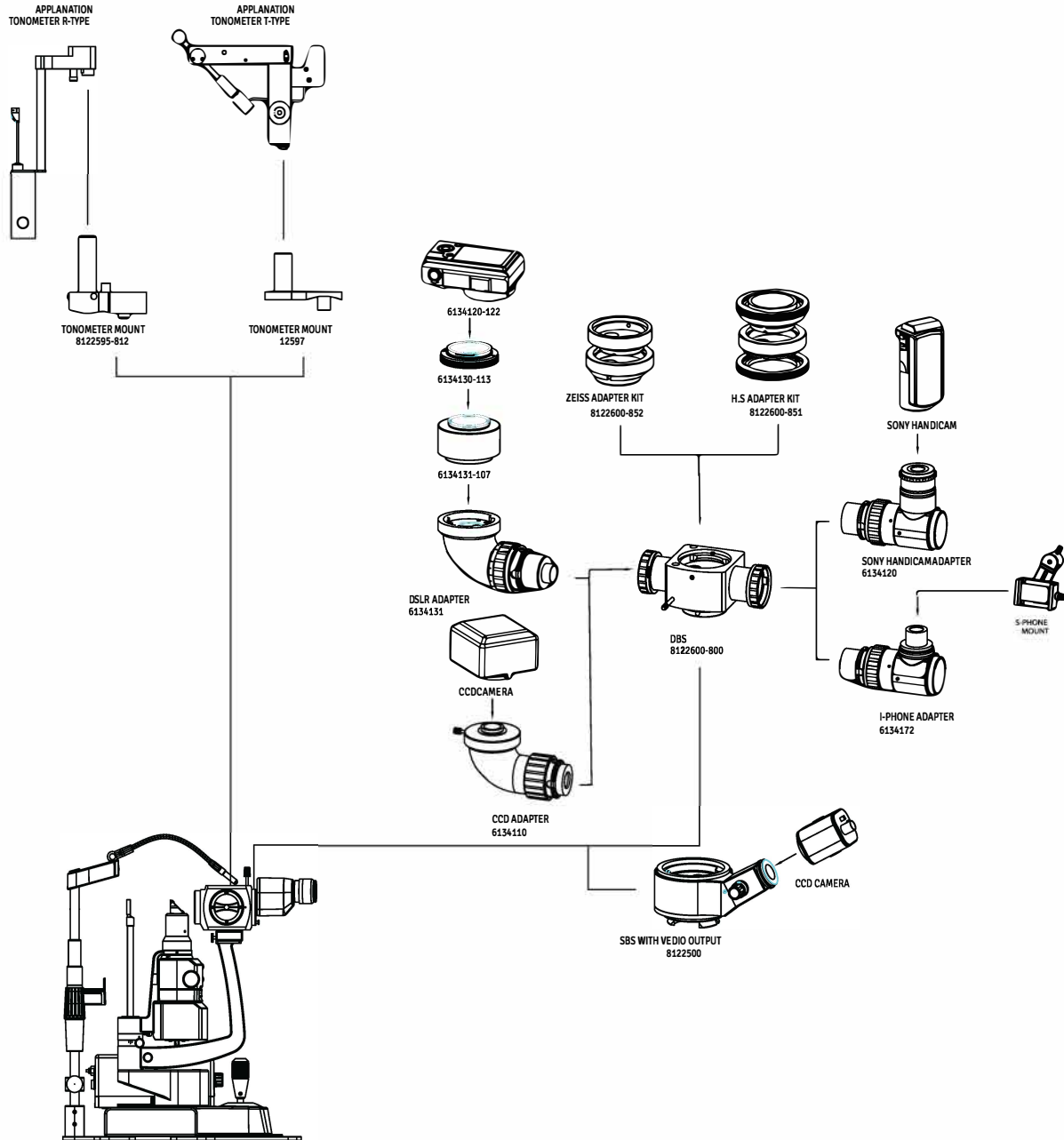
Chart of Common Errors

ISSUE	PROBABLE CAUSE	POSSIBLE SOLUTION
Lamp won't turn on.	Incorrect input power supplied to the eVO 300 Slit Lamp.	Check the outlet to ensure proper power is being supplied.
	Defective Power Cord.	Replace the Power Cord.
	Bulb may be blown out.	Replace Bulb.
	Defective Power Supply.	Replace the Power Supply.
Slit Lamp won't move.	Rubber stopper may be attached under the joystick.	Remove the rubber stopper.
	Base Lock Screw may be tightened.	Loosen the Base Lock Screw.
Rough base movement.	Rubber stopper may be attached under the joystick.	Remove the rubber stopper.
	Bearings may be damaged.	Replace the base.
	Shaft may be damaged.	Replace the base.
Fixation light does not light up.	Fixation Light Harness not plugged into the Power Supply Assembly.	Ensure the Fixation Light Harness is properly seated in the Power Supply Assembly.
	Defective Power Supply.	Replace the Power Supply.
Light too dim.	Incorrect wattage for bulb being used.	Replace with the proper Bulb.
	Bulb not installed properly.	Check bulb and ensure notch lines up with bulb housing.
Double slit visible in microscope.	Microscope not focused on focusing rod before use.	Install focusing rod and check to ensure microscope is focused on it.
	Bulb not installed properly.	Check bulb and ensure notch lines up with bulb housing.

The following is a checklist of items that need to be assessed in order to determine if the eVO 300 Slit Lamp requires Servicing.


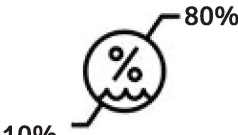
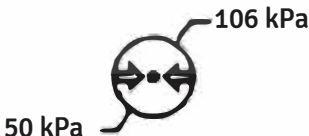
- Check the outside of the Slit Lamp for any damage or missing components.
- Inspect the power cord for damage.
- Test the lamp by turning the lamp on and turning the light all the way to its brightest setting, and all the way down to its lowest setting.
- Check to ensure all switches are functioning properly.
- Check the filters by cycling through all the options.
- Check the Slit Wheel by cycling through all the options.
- Check the base movement.

SYSTEM CONFIGURATION



SPECIFICATIONS

Catalog Number 8126200-795

Physical Dimensions	
Size:	Weight, unpacked: 23.0 lbs. (10.4 Kg)
Height: 19.8 in. (50.2 cm)	Weight, packed: 52 lbs (23.64 Kg)
Width: 10.5 in. (26.7 cm)	
Depth: 14.0 in. (35.6 cm)	
Electrical	
Voltage: 100-240 v AC	
Power consumption (max): 50-75VA	
Frequency: 50/60 Hz	
Fuses: T2.5AH250v	
Halogen Bulb: (6v, 20W)/ LED 5W.	
Operational Conditions	
Environmental	
The environmental conditions are as follows:	
Operating:	
Temperature: 10° C (50° F) to 35° C (95° F)	
Relative Humidity: 30% to 75%	
Atmospheric Pressure: 80 kPa (23.6 in. Hg) to 106 kPa (31.3 in. Hg)	
Transportation & Storage:	
Temperature -20° C (-4° F) to +70° C (158° F).	
Relative Humidity: 10% to 80% (non-condensing)	
Atmospheric Pressure: 50 kPa (14.8 in. Hg) to 106 kPa (31.3 in. Hg)	
Exposure to extreme temperature conditions indicated above must not exceed 15 weeks.	
Microscope	Galilean
Mag Change	5 Step Drum Rotation
Eyepiece	12.5X
Mag Ratio	6.5X, 10X, 16X, 25X, 40X
Field of View	37.0, 24.0, 15.0, 9.5, 6.0
PD Range	49 - 78 mm
Diopter Adjustment	+/- 5
Slit Illumination	6v 20W Halogen/ 5W LED
Slit Width	0 - 13.5 mm
Slit Length	0 - 13.5 mm
Slit Apertures	0.6, 5.8, 9, 13.5 mm and continuous 1.5 to 12.0 mm
Slit Rotation	0° - 180°
Slit Inclines	5, -10, -15, -20 Degrees
Filters	Red Free, Heat Absorbing, Cobalt Blue

SPECIFICATIONS (continued)

MOVEMENT RANGES

Longitudinal (In/Out)	100mm
Lateral (Left/Right)	107mm
Vertical (Up/Down)	30mm
Chin rest Range	80mm
Table Dimensions	18 1/2" x 12 10/16" (465 mm x 316 mm)

DISPOSAL

This product does not generate any environmentally hazardous residues. At the end of its product life, follow your local laws and ordinances regarding the proper disposal of this equipment.

SOFTWARE REVISION

There is no software installed in this unit.

GUIDANCE TABLES

<p align="center">Table 201 – Guidance and Manufacturer’s Declaration Electromagnetic Emissions All Equipment and Systems</p>		
<p align="center">Guidance and Manufacturer’s Declaration – Electromagnetic Emissions</p>		
<p align="center">The eVO 300 is intended for use in the electromagnetic environment specified below. The customer or user of the eVO 300 should ensure that it is used in such an environment.</p>		
Emissions Test	Compliance	Electromagnetic Environment - Guidance -
RF Emissions CISPR 11	Group 1 Class A	The eVO 300 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in near by electronic equipment.
Harmonics IEC 61000-3-2	Class A	The eVO 300 is suitable for use in all establishments, other than domestic, and those directly connected to the public low-voltage power network that supplies building used for domestic purposes.
Flicker IEC 61000-3-3	Complies	


GUIDANCE TABLES (continued)

Table 202 – Guidance and Manufacturer’s Declaration Electromagnetic Immunity All Equipment and Systems			
Guidance and Manufacturer’s Declaration – Electromagnetic Immunity			
The eVO 300 is intended for use in the electromagnetic environment specified below. The customer or user of the eVO 300 should ensure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
ESD IEC 61000-4-2	±6kv Contact ±8kv Air	±6kv Contact ±8kv Air	Floors should be wood, concrete or ceramic tile. If floors are synthetic, the R/H should be at least 30%.
EFT IEC 61000-4-4	±2kv Mains ±1kv I/Os	±2kv Mains ±1kv I/Os	Main power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1kv Differential ±2kv Common	±1kv Differential ±2kv Common	Main power quality should be that of a typical commercial or hospital environment.
voltage Dips/Dropout IEC 61000-4-11	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	Mains power quality should be that of a typical commercial or hospital environment. If the user of the eVO 300 requires continued operation during power mains interruptions, it is recommended that the eVO 300 be powered from an uninterruptible power supply or battery.
Power Frequency 50/60Hz Magnetic Field IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic field should be that of a typical commercial or hospital Environment.

GUIDANCE TABLES (continued)

<p align="center">Table 204 – Guidance and Manufacturer’s Declaration Electromagnetic Immunity Equipment and Systems that are NOT Life-supporting</p>			
<p align="center">Guidance and Manufacturer’s Declaration – Electromagnetic Immunity</p>			
<p>The eVO 300 is intended for use in the electromagnetic environment specified below. The customer or user of the eVO 300 should ensure that it is used in such an environment.</p>			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 vrms 150 kHz to 80 MHz	(v1) = 3 vrms	<p>Portable and mobile Rf communications equipment should be used no closer to any part of the eVO 300, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended Separation Distance:</p> $d=(3.5/v1)(\text{Sqrt } P)$ $d=(3.5/E1)(\text{Sqrt } P)$ 80 to 800 MHz $d=(7/E1)(\text{Sqrt } P)$ 800 MHz to 2.5 GHz
Radiated RF IEC 61000-4-3	80 MHz to 2.5 GHz @ 3V/m	(E1) = 3 v/m	
<p>Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p>Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p> <p>* Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and Tv broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. The measured field strength in the location in which the ME Equipment or ME System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the ME Equipment or ME System.</p> <p>* Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V1] V/m.</p>			

Interference may occur in the vicinity of equipment marked with the following symbol.



GUIDANCE TABLES (continued)

Table 206 – Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the eVO 300 for ME Equipment and ME Systems that are NOT Life-supporting.

Guidance and Manufacturer’s Declaration - Electromagnetic Immunity

Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the eVO 300

The eVO 300 is intended for use in the electromagnetic environment in which radiated RF disturbance are controlled. The customer or user of the eVO 300 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF Communications Equipment (transmitters) and the eVO 300 as recommended below, according to the maximum output power of the communications equipment.

Max Output Power of Transmitter (W)	Separation (m) 150kHz to 80 MHz $d=(3.5/\sqrt{f})(\text{Sqrt } P)$	Separation (m) 80 to 800 MHz $d=(3.5/E1)(\text{Sqrt } P)$	Separation (m) 800MHz to 2.5GHz $d=(7/E1)(\text{Sqrt } P)$
0.01	0.1166	0.1166	0.2333
0.1	0.3689	0.3689	0.7378
1	1.4666	1.4666	2.3333
10	3.6893	3.6893	7.3786
100	11.6666	11.6666	23.3333

For transmitters rated at a maximum output power not listed above, the recommended separation distance (d) in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (w) according to the transmitter manufacturer.
 Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.
 Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

WARRANTY

This product is warranted by Labo America Inc. against defective material and workmanship under normal use for a period of one year from the date of invoice to the original purchaser.(An authorized dealer shall not be considered an original purchaser). Under this warranty, Labo America Inc. sole obligation is to repair or replace the defective part or product at Labotech/Labomed discretion.

This warranty applies to new product and does not apply to a product that has been tampered with, altered in any way, misused, damaged by accident or negligence, or which has had the serial number removed, altered or effaced. Nor shall this warranty be extended to a product installed or operated in a manner not in accordance with the applicable LABOMED instruction manual, nor to a product which has been sold, serviced, installed or repaired other than by a Labo America Inc. factory or authorized LABOMED Dealer.

Lamps, bulbs, charts, cards and other expendable items are not covered by this warranty.

All claims under this warranty must be in writing and directed to the LABOMED factory, or authorized instrument dealer making the original sale and must be accompanied by a copy of the purchaser's invoice.

This warranty is in lieu of all other warranties implied or expressed. All implied warranties of merchantability or fitness for a particular use are hereby disclaimed. No representative or other person is authorized to make any other obligations for a LABOMED product. Labotech/Labomed shall not be liable for any special, incidental, or consequent damages for any negligence, breach of warranty, strict liability or any other damages resulting from or relating to design, manufacture, sale, use or handling of the product.

PATENT WARRANTY

If notified promptly on writing of any action brought against the purchaser based on a claim that the instrument infringes a U.S. Patent, Labo America Inc. will defend such action at its expense and will pay costs and damages awarded in any such action, provided that Labo America Inc. shall have sole control of the defense of any such action with information and assistance (at Labo America Inc. expense) for such defense, and of all negotiation for the settlement and compromise thereof.

PRODUCT CHANGES

Labo America Inc. reserve the right to make changes in design or to make additions to or improvements in its products without obligation to add such to product previously manufactured.

CLAIMS FOR SHORTAGES

We use extreme care in selection, checking, rechecking and packing to eliminate the possibility of error. If any shipping errors are discovered:

1. Carefully go through the packing material to be sure nothing was inadvertently overlooked when the unit was unpacked.
2. Call the dealer you purchased the product from and report the shortage. The material are packed at the factory and none should be missing if the box has never been opened.
3. Claims must be filed within 30 days of purchase.

CLAIMS FOR DAMAGES IN TRANSIT

Our shipping responsibility ceases with the safe delivery in good condition to the transportation company. Claims for loss or damage in transit should be made promptly and directly to the transportation company. If, upon delivery, the outside of the packing case shows evidence of rough handling or damage, the transportation company's agent should be requested to make a "Received in Bad Order" notation on the delivery receipt. If within 48 hours of delivery, concealed damage is noted upon unpacking the shipment and no exterior evidence of rough handling is apparent, the transportation company should be requested to make out a "Bad Order" report. This procedure is necessary in order for the dealer to maintain the right of recovery from the carrier.

Revision History

Rev. No.	Date of Release	DCR #	Change	App. By
1.4	April 4,2019	DCR/11/19	System Diagram added	S Bal
1.5	March 29,2023	ECR/07/2023	Electrical Specification Change	S Bal
1.6	Aug, 2024	DCR/11/2024	Regulatory Information	S Bal



Labo America Inc.
920 Auburn Court
Fremont, CA
94538
U.S.A.

Phone: 510-445-1257

Fax: 510-991-9862

Email: sales@laboamerica.com

www.laboamerica.com

EU|REP.

Labomed Europe
Essebaan 50
NL-2908 LK Capelle a/d IJssel
The Netherlands
Tel: +31 (0)10 4584222
Fax: +31 (0)10 4508251
E-mail: info@labomedeuropa.com



DQS Inc.
Labo America Inc.
Registered to ISO 9001:2015
BR Number 10005565



LABOTECH MICROSCOPES INDIA Pvt Ltd